Image Management Hub (IMH-20)
Quick Reference Guide

Starting and Ending a Case

1 Enter patient information

Select dropdowns and use virtual keyboard.




2 Record during the case

Start/Stop Recording (Options):

- Press scope button (if applicable).
- Press “REC” on IMH-20 touch screen.


Capture Image (Options):

- Press scope button (if applicable).
- Press “Capture” on IMH-20 touch screen.

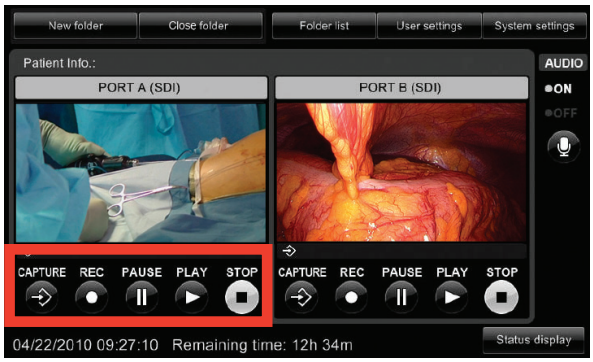


3 Ending the case

Select “Close Folder” then select “yes”.

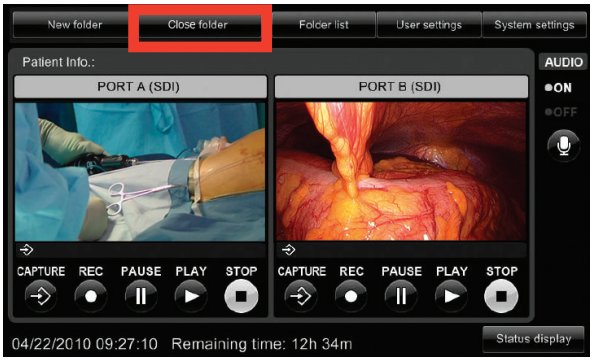

4 Starting the next case

Select “New Folder”.

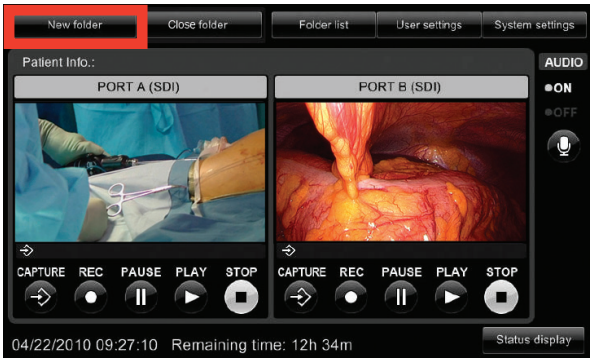


OLYMPUS®

Your Vision, Our Future


Image Management Hub (IMH-20)
Quick Reference Guide

Printing Images After a Case

1 Select “Folder List”

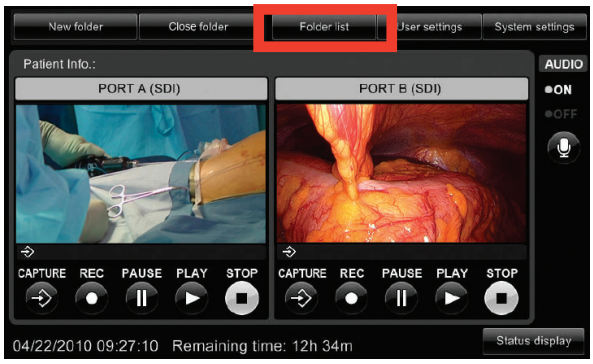


2 Select a file

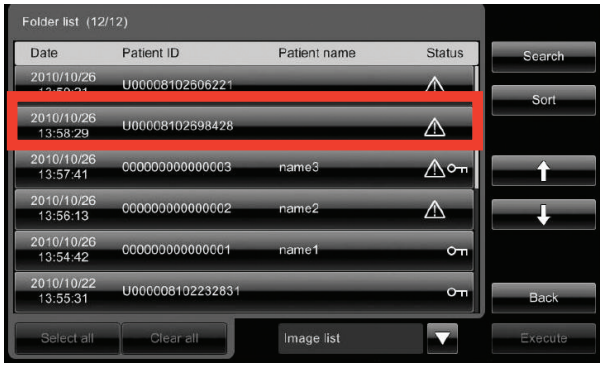


3 Select images

Highlight the images.

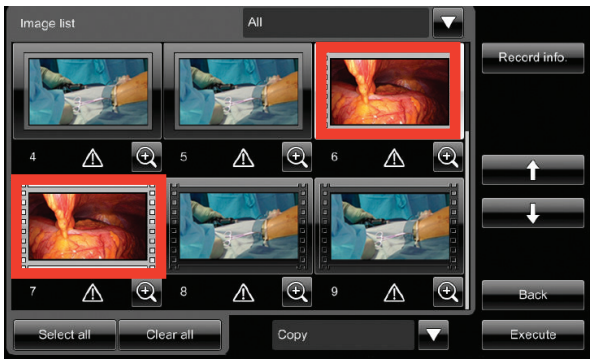


4 Select “Print”

Select “Print” from the dropdown menu.

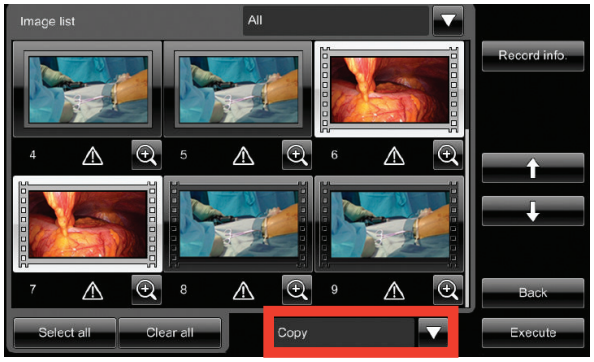

5 Select “Execute” to print

Select “Execute” to print the images.

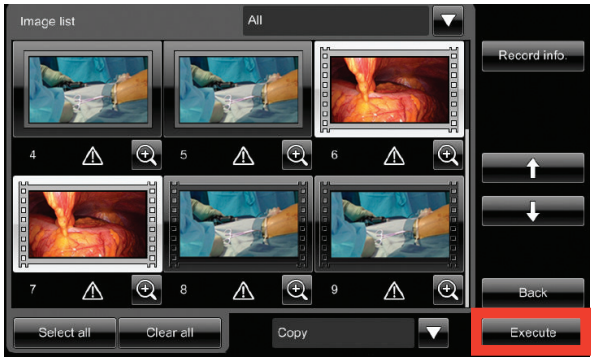


www.olympussurgical.com
©2012 Olympus America Inc. All rights reserved.
Printed in USA OAILAP0212QRG8678

OLYMPUS®
Your Vision, Our Future


Image Management Hub (IMH-20)
Quick Reference Guide

Copying to Flash Drive or External Hard Drive

1 Select “Folder List”

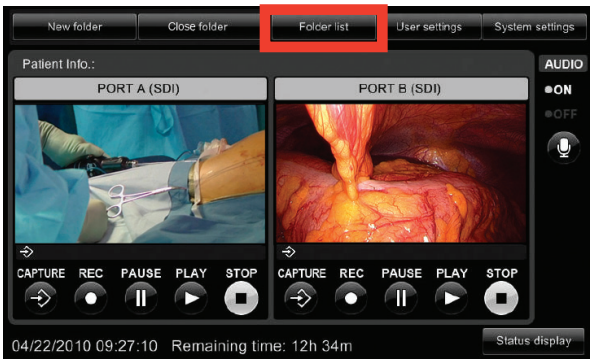


2 Select a file

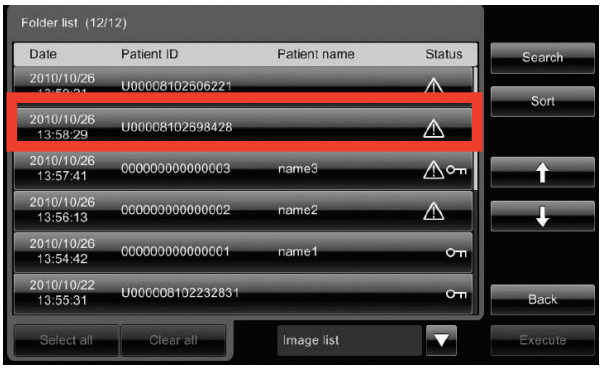


3 Select images and video

Highlight the images and video.

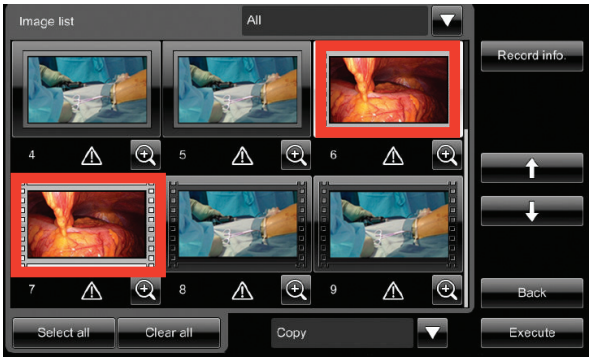

4 Select “Copy”


Select “copy” from the dropdown menu.

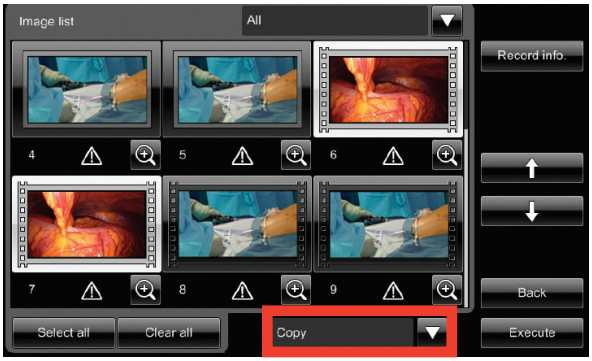

5 Insert a Flash Drive or External Hard Drive into “System 1”

6 Select “Execute” to copy

Select “Execute” to copy the images and video.

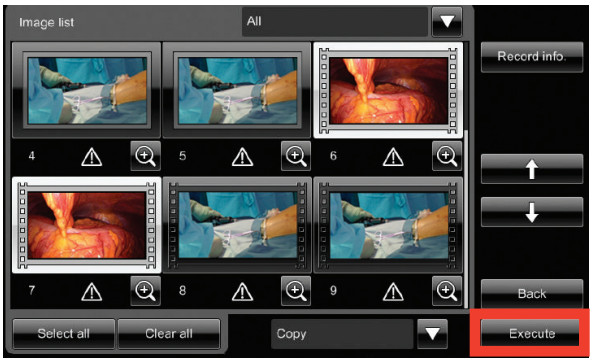


OLYMPUS®

Your Vision, Our Future


Image Management Hub (IMH-20)
Quick Reference Guide

Copying to DVD or Blu-Ray Disk

1 Select “Folder List”

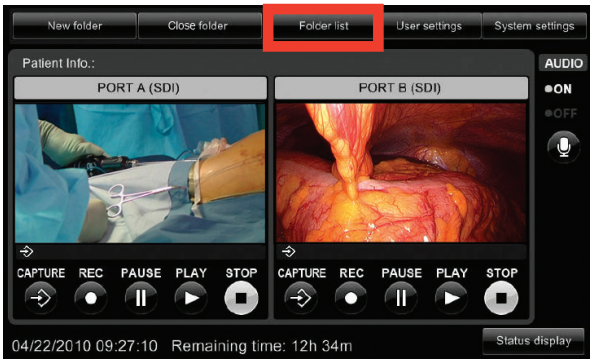


2 Select a file

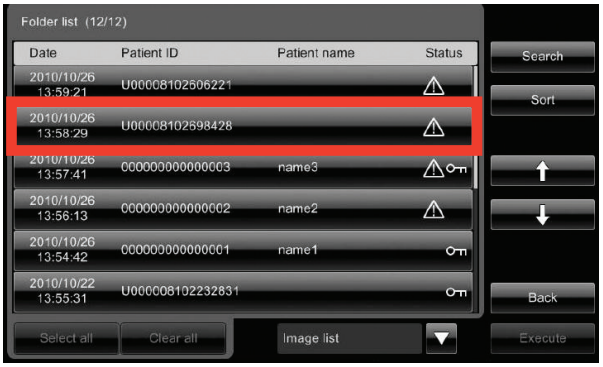


3 Select images and video

Highlight the images and video.

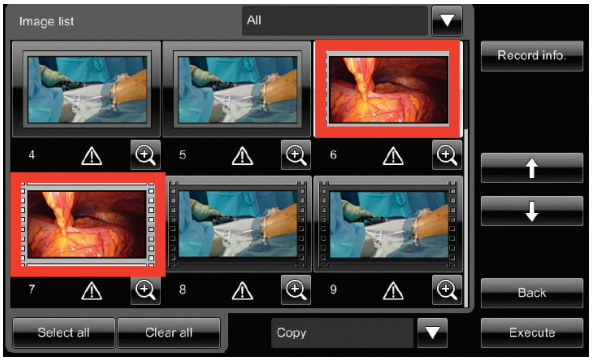


4 Select “Copy”

Select “copy” from the dropdown menu.

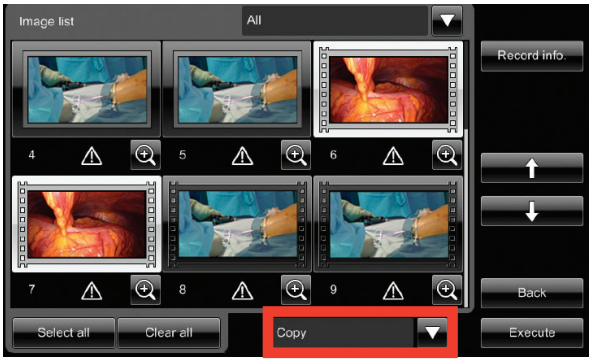

5 Select “Execute” to copy

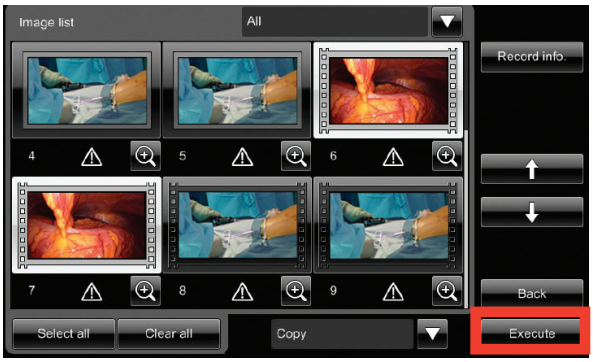

6 Insert a blank Blu-Ray disk or DVD, select “BD/DVD” from dropdown menu, then select “Execute”

Please wait a few seconds for the IMH-20 to recognize the disk.

www.olympussurgical.com
©2012 Olympus America Inc. All rights reserved.
Printed in USA OAILAP0212QRG8678

OLYMPUS®
Your Vision, Our Future