


## System Description

- Flat substrate, xyz stage, inkjet deposition system
- Low cost, user-fillable piezo-based inkjet print cartridges
- Built-in drop jetting observation system
- Fiducial camera for substrate alignment and measurement
- Variable jetting resolution and pattern creation PC-controlled with Graphical User Interface (GUI) application software
- Capable of jetting a wide range of fluids
- Heated vacuum platen
- Cartridge cleaning station
- Includes PC, monitor, and software



## Mechanical System


- Printable area
  - Substrate < 0.5 mm thickness: 210 mm x 315 mm (8.27 in x 12.4 in)
  - Substrate 0.5 - 25 mm thickness: 210 mm x 260 mm (8.27 in x 10.2 in)
- Repeatability:  $\pm 25 \mu\text{m}$  ( $\pm 0.001$  in)
- Substrate holder
  - Vacuum platen
  - Temperature adjustable; ambient to 60° C
- System Footprint: 673 mm x 584 mm x 419 mm (26 in x 23 in x 16 in)
- Weight approximately 43 kg (95 lbs)
- Power 100-120/200-240 VAC 50/60 Hz 375 W maximum
- Operating range 15-40° C at 5-80% RH non-condensing
- Altitude up to 2000 m
- Safety and EMC compliance
  - Safety: NRTL Certified to EN 61010-1, UL 61010-1, CSA 22.2 No. 61010-1
  - EMC: EN61326-1 Class A, FCC Part 15 Class A

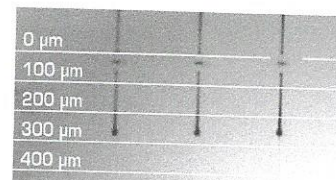

Drop Watcher View

## Fiducial Camera

- Allows substrate alignment using reference marks
- Allows positioning a print origin or reference point to match substrate placement
- Provides measurement of features and locations
- Provides inspection and image capture of printed pattern or drops
- Provides cartridge alignment when using multiple cartridges
- Allows matching drop placement to previously patterned substrate

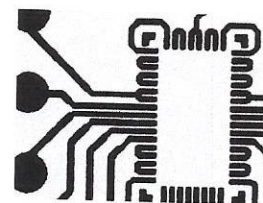

Fiducial Camera View

## Cartridge

- Type: Piezo-driven jetting device with integrated reservoir and heater
- Usable Ink Capacity: Up to 1.5 ml (user-fillable)
- Materials Compatibility: Many water-based, solvent, acidic or basic fluids
- Number of Nozzles: 16 nozzles, 254  $\mu\text{m}$  spacing, single row
- Drop Volume: 1 (DMC-11601) and 10 (DMC-11610) picoliter nominal



## Control PC and Application Software

- Pre-loaded patterned templates
- Pattern preview
- Editors: Pattern, piezo drive waveform, cleaning cycle, substrate setting
- Bitmap (1 bit) files accepted
- DXF, Gerber, GDSII and OASIS file conversion to Bitmap


## Replaceable Items

- Print cartridge with one-time user-fillable reservoir
- Cleaning station nozzle blotting pad
- Drop watcher fluid absorbing pad

**Dimatix**

TM



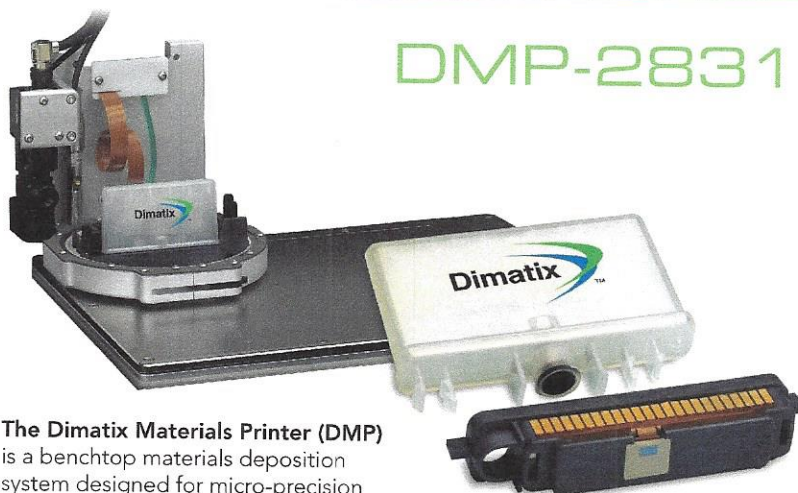

## DMP-2831

**The Dimatix Materials Printer (DMP)** is a benchtop materials deposition system designed for micro-precision jetting a variety of functional fluids onto virtually any surface, including plastic, glass, ceramics, and silicon, as well as flexible substrates from membranes, gels, and thin films to paper products. A complete turnkey system, the DMP facilitates developing and testing manufacturing processes and product prototypes. It also can be used for prototyping of products from flexible circuits, RFID tags and displays to DNA arrays, and wearable electronics. By employing inexpensive exchangeable cartridges that researchers can fill with their own fluid materials, the DMP system minimizes waste of expensive fluid materials, thereby eliminating the cost and complexity associated with traditional product development and prototyping.

**The Dimatix Materials Cartridge** is a cartridge-based inkjet printhead used with the DMP and available in 1 pl and 10 pl drop volumes. Based on FUJIFILM Dimatix's proprietary Silicon MEMS technology, the 16-jet Dimatix Materials Cartridge is designed for high-resolution, non-contact jetting of functional fluids in a broad range of applications. The industry-first 1 pl cartridge can deposit features as small as 20  $\mu\text{m}$  (20 millionths of a meter) to fabricate products such as organic thin-film transistors (TFTs) and printed circuits. In biotechnology, the Dimatix Materials Cartridge allows researchers to closely pack large numbers of elements in DNA arrays, to permit more accurate and efficient analyses.

# FUJIFILM



**Corporate Office:**  
FUJIFILM Dimatix, Inc.  
2250 Martin Avenue  
Santa Clara, CA 95050  
USA

Tel: (408) 565-9150  
Fax: (408) 565-9151  
Email: info@dimatix.com

**New Hampshire Facility:**  
FUJIFILM Dimatix, Inc.  
109 Etna Road  
Lebanon, NH 03766  
USA

Tel: (603) 443-5300  
Fax: (603) 448-9870  
Email: info@dimatix.com

**Japan Office:**  
FUJIFILM Global Graphic Systems Co., Ltd.  
2-26-30 Nishiazabu  
Minato-ku, Tokyo 106-0031  
Japan  
Advanced Marking Strategy Division  
Phone: +81 3 6419 0530  
Fax: +81 3 6419 9840  
E-mail: dmp\_ffgs@ffgs.fujifilm.co.jp

**Europe Office:**  
Tel: +44 7739 863 505  
Fax: +44 870 167 4328  
Email: euro@dimatix.com

**Korea Office:**  
Email: mdkorea@dimatix.com

**Taiwan Office:**  
Email: mdtaiwan@dimatix.com

**China Office:**  
FUJIFILM Dimatix China Service Center  
Building 30, 1000 Jinhai Road  
Pudong New Area, Shanghai  
China 201206  
Email: china@dimatix.com

**Singapore Office:**  
Email: mdsingapore@dimatix.com