



Macofar

CD 20 CD 5 CAPSULE FILLING



- ◆ Powder, pellets, tablets
- ◆ Proven reliability
- ◆ Powder containment
- ◆ Fastest format change-over




The Macofar compact, intermittent-motion capsule filling and closing machines, CD5 and CD20 are ideal for laboratory, clinical-trial and small production (mechanical speed of 6,000 and 20,000 capsules/hour).

The main features of this machine are:

- ◆ outstanding filling accuracy of vertical-movement vacuum dosators;
- ◆ capability to fill a wide range of powders and/or pellets and/or tablets;
- ◆ limited format parts for fastest format change-over;
- ◆ proven field reliability, easy to use, and economic to run;
- ◆ full technical documentation for FDA validation;
- ◆ CE and GMP conform.

CD 20 CD 5

CAPSULE FILLING


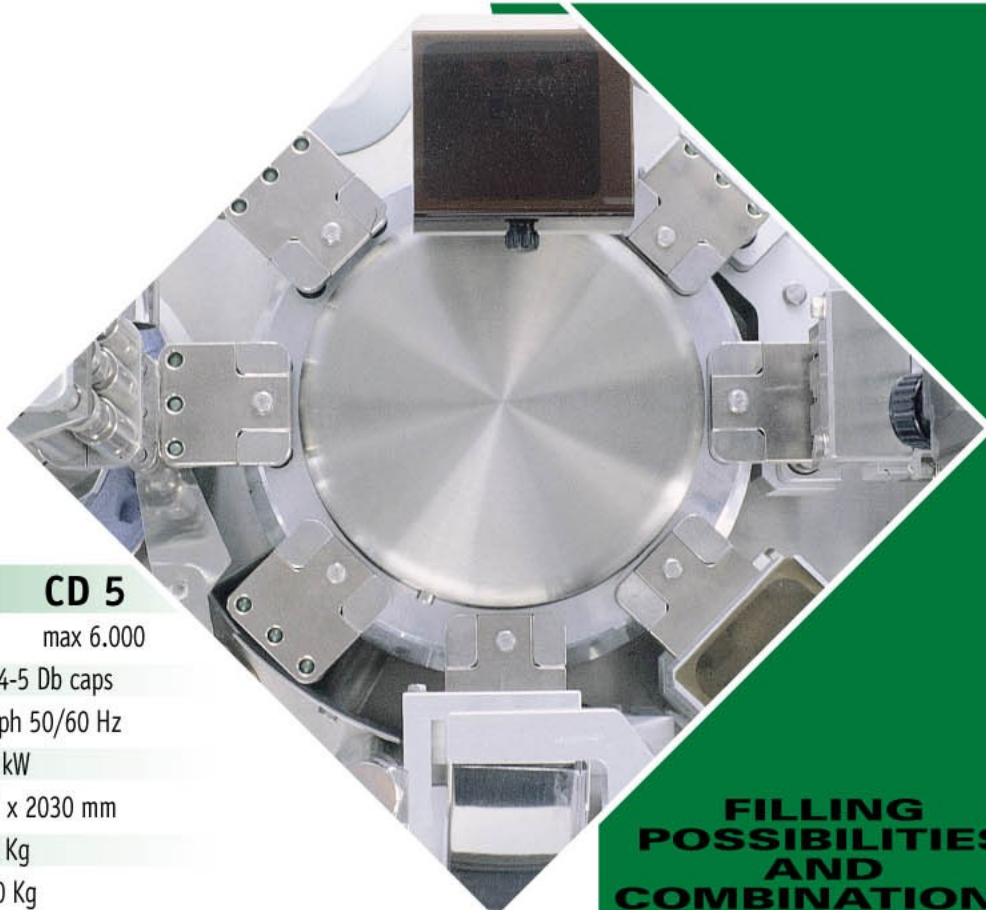


Technical data

CD 20 CD 5

Output capsules/hour	max 20.000	max 6.000
Capsule size	00-0-1-2-3-4-5 Db caps	
Standard Voltage	220/380 V-3ph 50/60 Hz	
Rated Power	2,5 kW	
Packing	1970 x 1670 x 2030 mm	
Net weight	800 Kg	
Gross weight	1150 Kg	
Machine compressed air	100 NL/1' a 6 bar	50 NL/1' a 6 bar
Machine + C.W. compressed air	300 NL/1' a 6 bar	200 NL/1' a 6 bar

Main Options
 CS60 capsule sorter
 CW statistic weigher
 W100% weight control

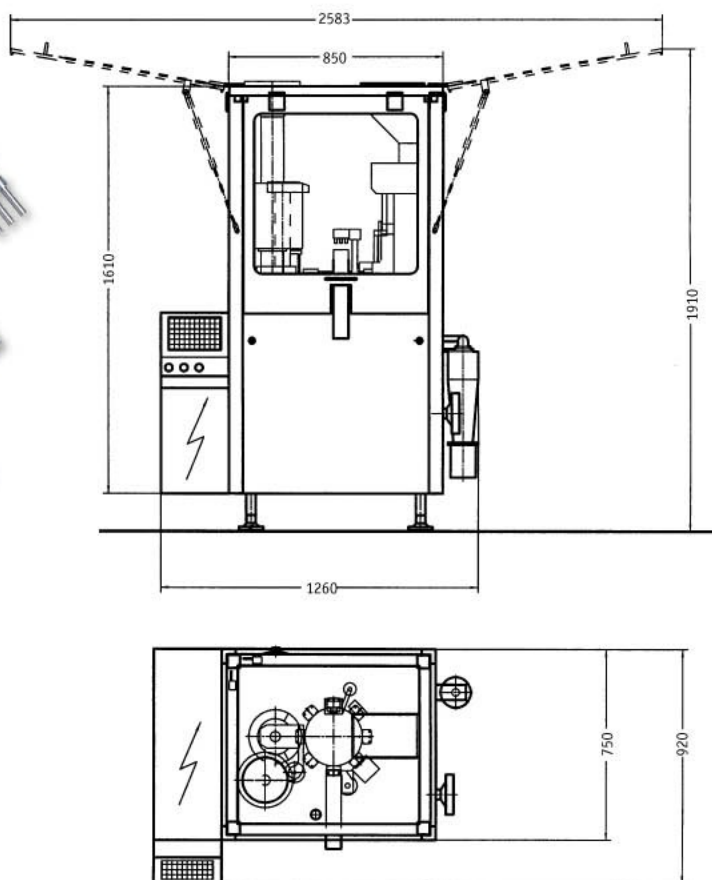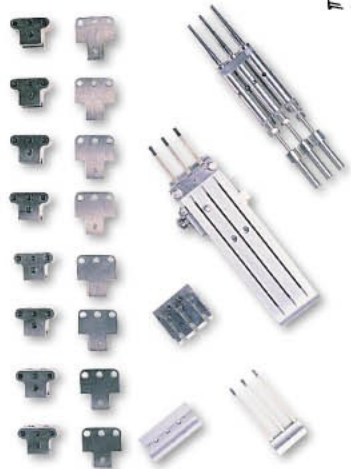

FILLING POSSIBILITIES AND COMBINATIONS

